



Mechanical Circulatory Support LVAD, RVAD, and Bi-VAD



History

- SAMPLE
- Bridge to transplant
- Destination therapy
- Estimated downtime
- LVAD, RVAD, Bi-Vad, TAH
- DNR, MOST, or Living Will
- Contact with LVAD coordinator

Signs and Symptoms

- Unconsciousness
- Pulseless
- Apneic
- Poor capillary refill / skin color
- AMS or decreased mental status
- No electrical activity on ECG
- No heart tones on auscultation

Differential

- See Reversible Causes below
- Infection/Sepsis
- Hypovolemia
- Cardiac arrest
- Hemorrhage

Contact VAD coordinator:

- As quickly as possible for troubleshooting and treatment advice, but do not delay emergency treatment
- Follow patient specific emergency plan if present

Rapid assessment
Check for signs of life
Assess for adequate perfusion

**Criteria for Death /
No Resuscitation**
Review DNR / MOST Form

NO

Decomposition
Rigor mortis
Dependent lividity
Blunt force trauma
Injury incompatible with life
Extended downtime with asystole

Do not begin resuscitation

Follow
Deceased Subjects
Policy

**Unresponsive and
Not breathing normally**

Assess LVAD function
Look and listen for alarms
LVAD Alarming?

NO

Place stethoscope over heart

Humming sound present?

NO

Airway Protocol(s) AR 1, 2, 3
if indicated

Respiratory Distress Protocol AR 4
if indicated

Assume VAD failure
Initiate age appropriate ventilation rate

Go to Page 2

**Responsive or Unresponsive and
Breathing normally**

Airway Protocol(s) AR 1, 2, 3
if indicated

Respiratory Distress Protocol AR 4
if indicated

Assess LVAD function
Look and listen for alarms
LVAD Alarming?

NO

Place stethoscope over heart

Humming sound present?

NO

Assume VAD failure
Initiate age appropriate ventilation rate

Go to Page 2

Go to Page 2

Go to Page 3



Mechanical Circulatory Support LVAD, RVAD, and Bi-VAD



History


- SAMPLE
- Bridge to transplant
- Destination therapy
- Estimated downtime
- LVAD, RVAD, Bi-Vad, TAH
- DNR, MOST, or Living Will
- Contact with LVAD coordinator

Signs and Symptoms

- Unconsciousness
- Pulseless
- Apneic
- Poor capillary refill / skin color
- AMS or decreased mental status
- No electrical activity on ECG
- No heart tones on auscultation


Differential

- See Reversible Causes below
- Infection/Sepsis
- Hypovolemia
- Cardiac arrest
- Hemorrhage

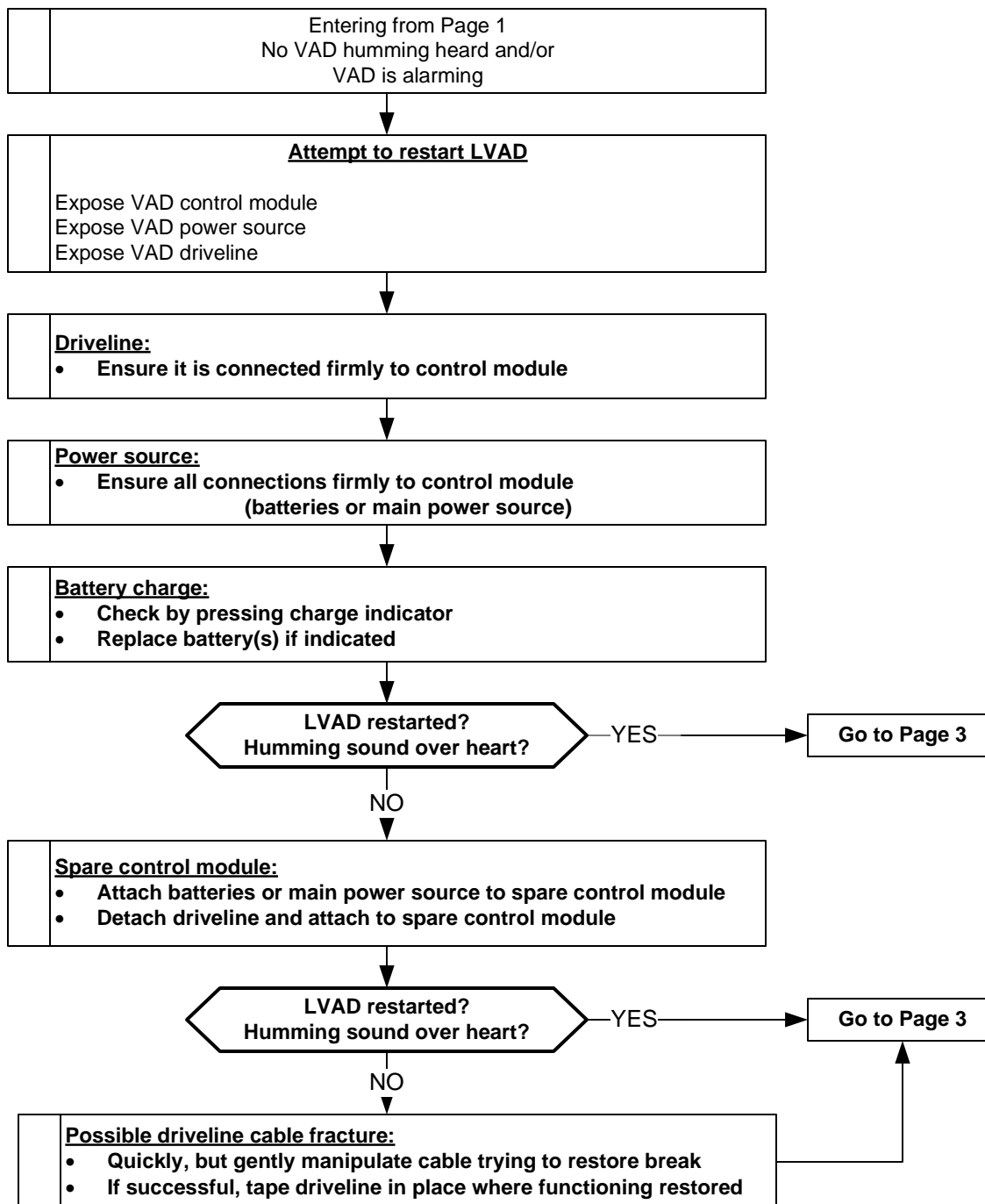







Mechanical Circulatory Support LVAD, RVAD, and Bi-VAD

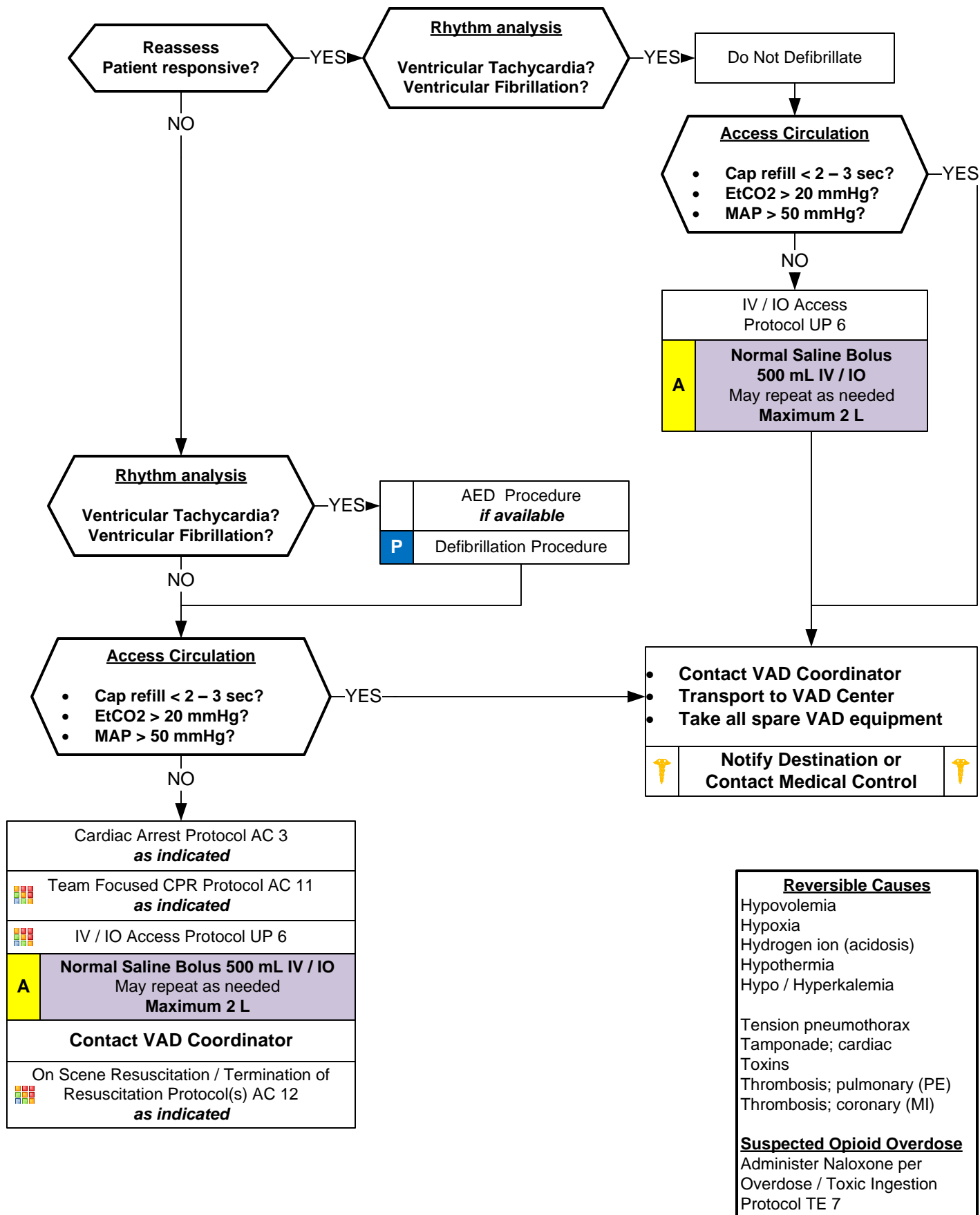

HeartMate II

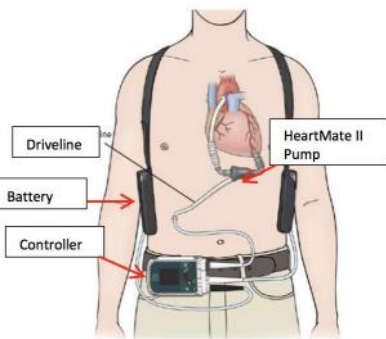

HeartMate 3

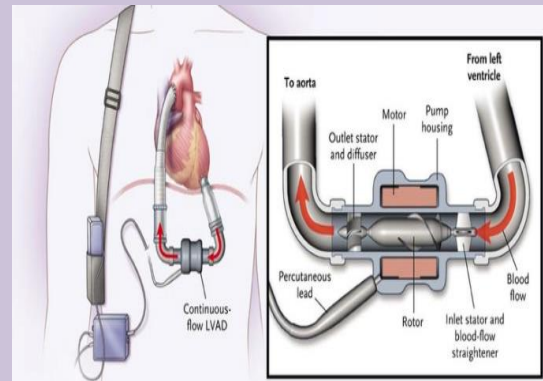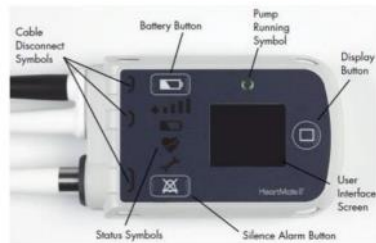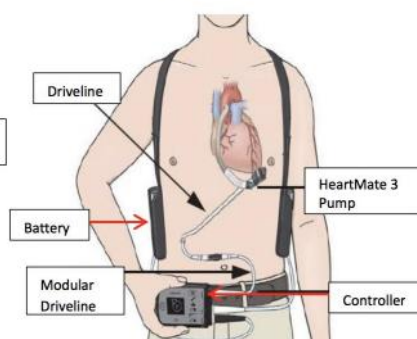

Pearls

- **Recommended exam:** Mental status, skin color, capillary refill, peripheral pulses, blood pressure.
- **Assessment of blood flow and perfusion status:**
 - Optimal BP attained by manual BP and Doppler.
 - Automated BP devices can measure a BP in about 50% of attempts and is not reliable to assess perfusion
 - A MAP of ≥ 60 mmHg is adequate for most LVAD patients.
 - Skin color, skin temperature, capillary refill
- **Mechanical Circulatory Support devices:**
 - LVAD – Left Ventricular Assist Device
 - RVAD – Right Ventricular Assist Device
 - BiVAD – Biventricular Ventricular Assist Device
 - TAH – Total Artificial Heart
- **Reasons for use:**
 - Bridge therapy – patients awaiting transplant or anticipated recovery.
 - Destination therapy – advanced heart failure, not candidate for transplant, and will live rest of life with device.
- **Pump type and assessing pulses:**
 - Pulsatile flow pumps – older units, not commonly in use now, but generate blood flow with a pulsatile flow and patient will have a palpable pulse.
 - Continuous flow pumps – majority of pumps now used and create blood flow in a continuous stream, no pulsatile flow, so patient will not have a palpable pulse.
 - Most devices are implanted inside the chest and have an internal pump, a driveline connected from the pump to the controller unit, and a power source consisting of batteries and electrical cord for receptacles.
- **Common complications:**
 - Disconnection of power supply, either battery disconnect, or electrical cord to receptacle disconnection.
 - Driveline failure or disconnection from controller unit.
 - Controller failure
 - Blood clot formation, acute stroke, and bleeding (mucosal and gastrointestinal most common sites)
 - Infection
- **Abnormal heart rhythm:**
 - Pseudo-PEA: Normal cardiac electrical activity in a patient who is alert and well perfused with no palpable pulse.
 - Tachyarrhythmias are usually well tolerated.
 - Ventricular fibrillation and tachycardia are treated in the usual fashion.
- **End Tidal CO₂ (EtCO₂)**
 - If EtCO₂ is < 10 mmHg, improve chest compressions. Goal is ≥ 20 mmHg.
 - If EtCO₂ spikes, typically > 40 mmHg, consider Return of Spontaneous Circulation (ROSC)
- **Transcutaneous Pacing:**
 - Pacing is NOT effective in cardiac arrest and pacing in cardiac arrest does NOT increase chance of survival